



Safety at every step

CareMat

The pressure-sensitive sensor mat for healthcare.



Making daily life easier in care

In care facilities, the CareMat plays a significant role in maintaining residents' quality of life and also provides valuable support to care staff.

The CareMat is simply placed in **front of the bed or the room door** and connected via **cable or wireless** to either the facility's call system or a wireless receiver plugged into an outlet.

1



People with dementia tend to be restless and disoriented. They often try to get up without assistance.

2



In this case, the CareMat provides reliable protection. Stepping on the CareMat triggers either the **call signal from the home's nurse call system** or a **visual or audible signal from the outlet receiver**. Caregivers can quickly provide assistance.

3



The CareMat makes the work of nursing staff easier and increases patient safety. It helps prevent falls and stops patients from leaving the room unsupervised.

4



CareMat is your trusted partner in supporting caregivers and residents.



Use of CareMat

Caring for people with dementia or those at risk of falling or wandering poses a major challenge for caregivers. The focus is on avoiding unnecessary restrictions on the mobility of those in need of care.

The pressure-sensitive CareMat sensor mat provides valuable support in this regard.

Areas of application

- Retirement & nursing homes
- Hospitals & clinics
- Psychiatric facilities
- Assisted living communities
- In-home care

Easy installation

- Connection to the building's intercom system
- Connection to a wall-mounted receiver
- Placement in front of the bed or the room door

Features

- Reliable alarm
- Portable
- Non-slip surface
- Hygienic material
- Durable construction

Benefits for caregivers

- Prevention of falls
- Prompt assistance
- Reduction in nightly rounds
- Low follow-up costs
- Easy to use

Benefits for patients


- Unrestricted freedom of movement
- Sense of security
- Improved quality of life

Technical specifications

CareMat can be placed in front of the care bed or the door and is available in wired or wireless version.

New: The wireless transmitter module is now even easier to replace!

Simply pull the magnetic wireless transmitter module out of its mount, slide the wireless circuit board to the side, and replace the battery without using screws.

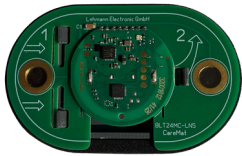

CareMat A



CareMat B



CareMat C

Dimensions

CareMat A - 1,100 x 700 x 9 mm, 8.2 kg

CareMat B - 700 x 400 x 9 mm, 3.2 kg

CareMat C - 1,100 x 700 x 9 mm, 6.5 kg

Material

Polyurethane

Surface

Nubbed texture

Color

Careyellow, Caregrey

Contact output

Cable: Active or idle contact, resistive circuit

Wireless: Manufacturer-specific wireless protocol

Connection

Assembled by your supplier to match your call system.

Your CareMat Distribution Partner

meta-care is your europe-wide distributor of medical devices and supplies for the nursing and healthcare sectors.

We offer a comprehensive range of products for a wide variety of care situations and are your expert partner for every equipment related to nursing beds. With meta-care, you can count on quality, safety, and practical innovations.

meta-care GmbH

T: +43 5 1721 0

E: office@meta-care.eu

W: www.meta-care.eu

Akazienhofstraße 20

9020 Klagenfurt am Wörthersee

Austria